



WHOLE STOMIZATION ACCEPT






Nowadays, in surgeries such as laparoscopic surgery, developing gauze towel or gauze swab with developing function are widely used.

The main function of the gauze towel or gauze swab is to absorb fluid and protect the organ tissue to avoid secondary damage.

At present, there is a solution for the use of developing cotton patties in neurosurgery.


Surgical patties are cotton sheets, which have better performance than gauze, and are more capable of applying to such dangerous operations.

## Feature

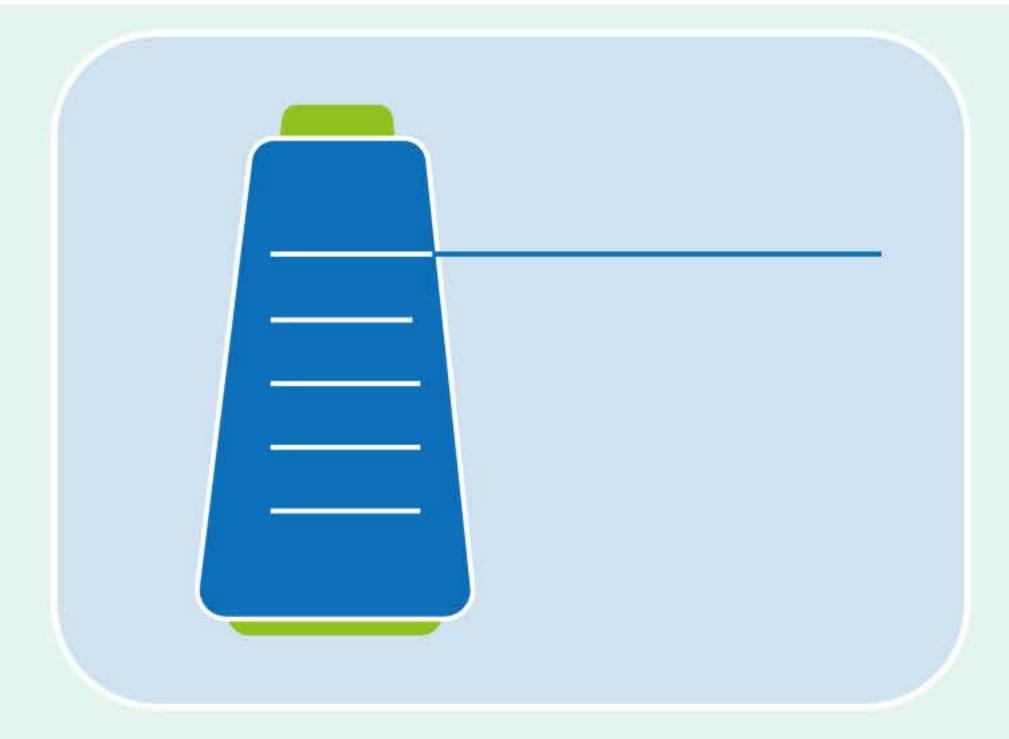


PVC/TPU Detectable Thread under X-ray scan

- · 0.3-0.5mm thickness thread
- · Detectable but keep a fabric feature
- · Special formula to make it firm and strong than usual

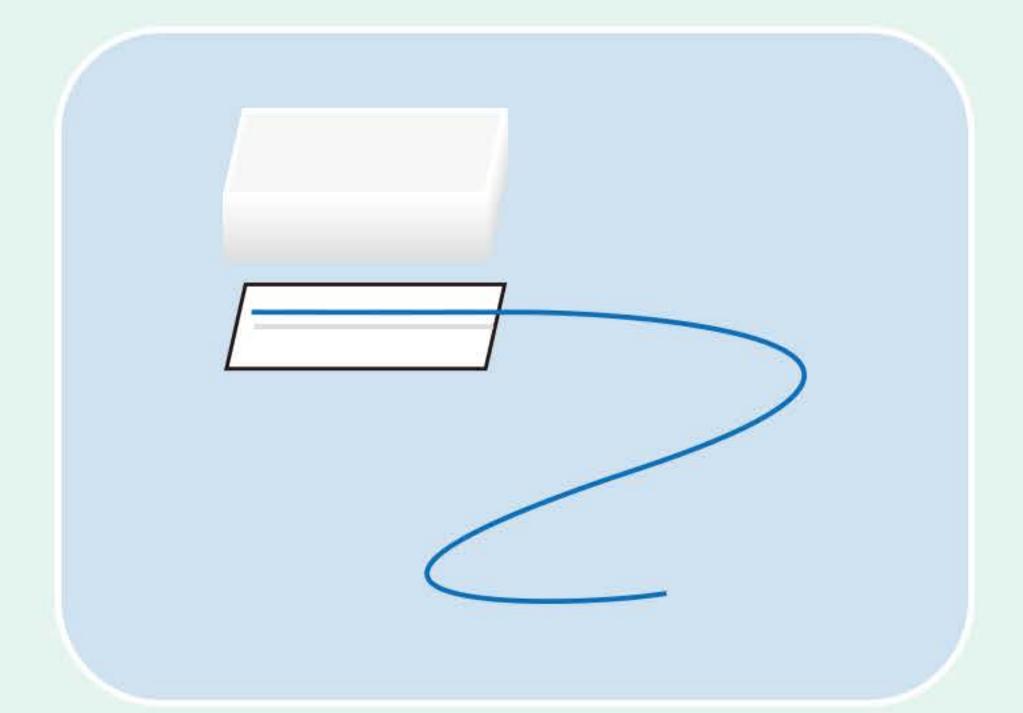

Integrated Manufacturing process

- Whole thread detectable
- · Ultrasonic bonding to avoid falling and sheding
- No toxic and impurities

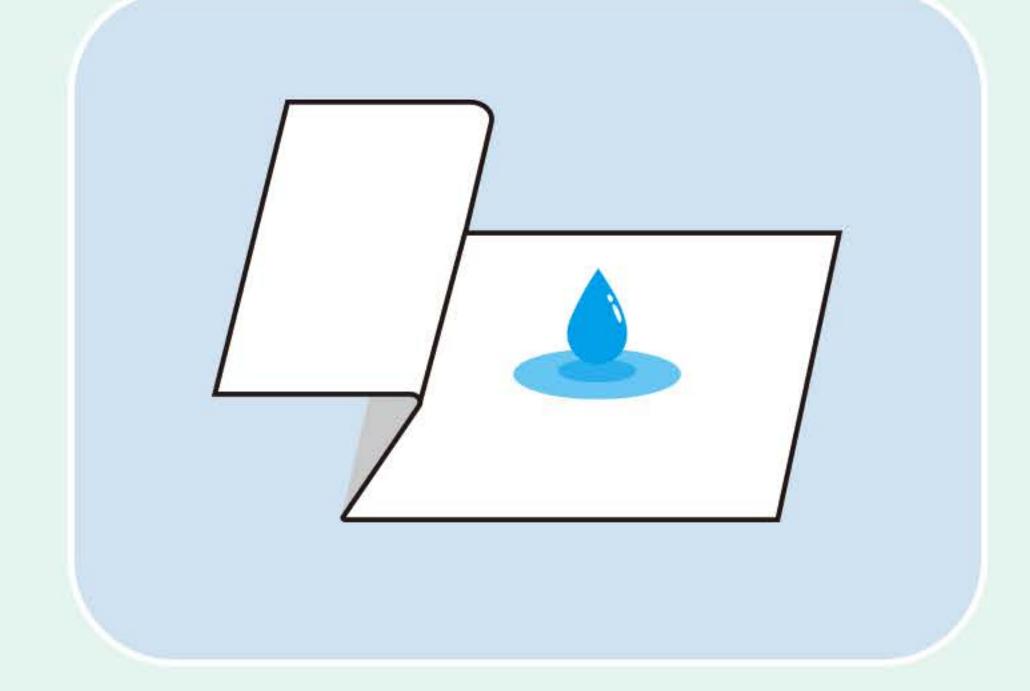

Spunlace 100% cellulose cotton

- Spunlace process can make fiber tighter, and cause no pollution or damage. The cotton has no adhesive, no pilling, no shedding, and no other impurities.
- · Spunlace cellulose cotton features moisture absorption, softness, high strength, good appearance and touch
- · 3-4 times fluid absorbency capacity based on itself.

## NANJING SENGONG BIOTECHNOLOGY CO.,LTD

B8-201,#9 Kechuang Avenue, Jiangbei New Area, Nanjing City, Jiangsu, 210000, China



www.purnote.com



Grel\_xu@purnote.com





