

TPM3-NTRK1 [G667C]/BaF3

CBP73210

Contents

I. Introduction.....	1
II. Background.....	1
III. Representative Data.....	2
1. WB of TPM3-NTRK1[G667C]/BaF3	2
2. Sanger of TPM3-NTRK1[G667C]/BaF3	2
3. Anti-proliferation assay	3
IV. Handling Procedure for Flask Cultures.....	3
V. Subculturing Procedure.....	4
VI. Cryopreservation Procedure.....	4
VII. Performing the Anti-proliferation Assay.....	5
1. Materials to Be Supplied by the User.....	5
2. Protocol for the Cell Viability Assay.....	5
VIII. References.....	6



TPM3-NTRK1 [G667C]/BaF3

CBP73210

I. Introduction

Cell Line Name: TPM3-NTRK1 [G667C]/BaF3

Host Cell: Ba/F3

Stability: 16 passages (in-house test, that not means the cell line will be instable beyond the passages we tested.)

Application: Anti-proliferation assay and PD assay

Freeze Medium: 90% FBS+10% DMSO

Complete Culture Medium: RPMI-1640+10%FBS+2ug/ml puromycin

Mycoplasma Status: Negative

II. Background

NTRK1 rearrangements involve the kinase domain of the NTRK1 gene which includes 17 exons located on chromosome 1q21-22 (Alberti et al. 2003). The frequency of NTRK1 fusions in patients with adenocarcinoma histology is 3.3% of cases (3 out of 91 patients; Vaishnavi et al. 2013).

a In preclinical studies, cell lines expressing MPRIP-NTRK1 or CD74-NTRK1 were inhibited by ARRY-470, a TRKA/B/C inhibitor (Vaishnavi et al. 2013).

b In preclinical studies, cell lines expressing MPRIP-NTRK1 or CD74-NTRK1 were inhibited by lestaurtinib, a FLT3/TRKA inhibitor (Vaishnavi et al. 2013).


c A patient with non-small cell lung cancer harboring an MPRIP-NTRK1 fusion mutation was treated with crizotinib, which has activity against



TRKA in addition to ALK, MET and ROS1 (Vaishnavi et al. 2013). The patient demonstrated a small radiographic decrease in tumor size but experienced disease progression after three months (Vaishnavi et al. 2013)

III. Representative Data

1. WB of TPM3-NTRK1-G667C/BaF3 expression

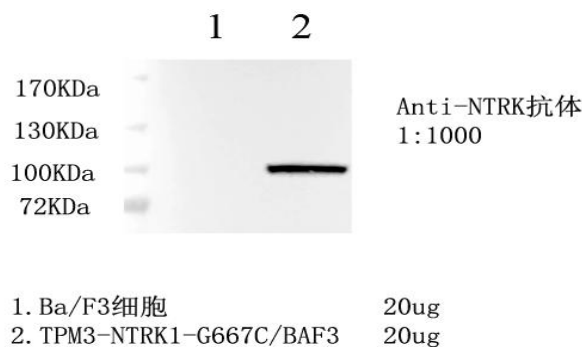

2. Sanger of TPM3-NTRK1 [G667C]/BaF3

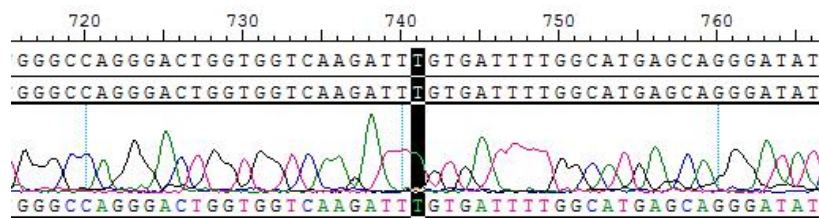


Figure 2. SNP of TPM3-NTRK1-G667C/BaF3

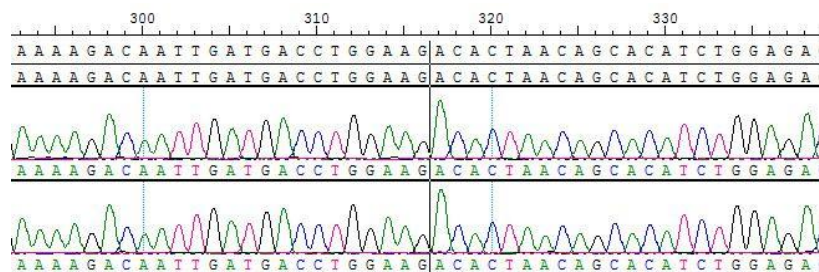


Figure 3. Sanger Sequencing of TPM3-NTRK1-G667C/BaF3



3. Anti-proliferation assay

CTG Proliferation Assay of BaF3 TPM3-NTRK1 G667C Cells (C1)

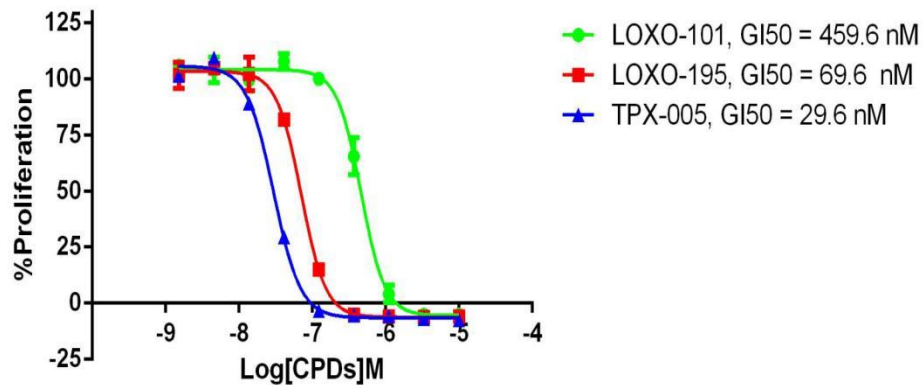


Figure 4. CTG Proliferation Assay of BaF3 TPM3-NTRK1 G667C Cells (C1).

IV. Handling Procedure for Flask Cultures

The flask was seeded with cells grown and completely filled with medium at Cobioer.

1. Upon receipt visually examine the culture for macroscopic evidence of any microbial contamination.

Using an inverted microscope (preferably equipped with phasecontrast optics), carefully check for any evidence of microbial contamination.

2. Incubate the flask in an upright position for several hours at 37°C.

After the temperature has equilibrated, aseptically remove the entire contents of the flask and centrifuge at 1000 rpm for 5 minutes. Remove shipping medium and save for reuse. Resuspend the cell pellet in 10 ml of this medium.

3. From this cell suspension remove a sample for a cell count and viability. Adjust the cell density of the suspension to 3×10^5 viable cells/ml



in the shipping medium.

4. Incubate the culture, horizontally at 37°C in a 5% CO₂ in air atmosphere. Maintain the cell density of the culture as suggested under the subculture procedure.

V. Subculturing Procedure

Cultures can be maintained by the addition of fresh medium or replacement of medium. Alternatively, cultures can be established by centrifugation with subsequent resuspension at 3*10⁵ viable cells/mL. Do not allow the cell density to exceed 2*10⁶ cells/mL.

Interval: Maintain cultures at a cell concentration between 3* 10⁵ and 2* 10⁶ viable cells/mL.

Medium Renewal: Add fresh medium every 2 to 3 days (depending on cell density)

VI. Cryopreservation Procedure

1. Transfer cell suspension(4-8*10⁶ cells) to a centrifuge tube and spin at approximately 1000rpm for 5 minutes.
2. Discard supernatant and resuspend cells in cryopreservation medium(90% FBS+10%DMSO).
3. Transfer the cells into Freezing Tube, 1ml/vial.
4. Place Freezing Tube in a Frosty container that is kept at room temperature and has sufficient isopropanol.
5. Place the Frosty container into the -80°C freezer overnight.
6. On next day, transfer the vials into liquid nitrogen.



VII. Performing the Anti-proliferation Assay

1. Materials to Be Supplied by the User

- Compounds to be tested
- RPMI-1640 (Gibco #C11875500CP)
- FBS (Gibco #10099-141C)
- DMSO
- CellTiter-Glo® Luminescent Cell Viability Assay (Promega, Cat.No.: G7573)
- 96 Well Assay Plate (White Plate, Clear Bottom with Lid Tissue Culture Treated Polystyrene 1/Pack, Corning #3610)
- T25 Flask
- 96 Well Storage Microplate(Corning#3357)
- Synergy H1 Hybrid Multi-Mode Reader(Biotek)

2. Protocol for the Anti-proliferation Assay

- 1) Take logarithmic growth cells, centrifuge and discard the culture supernatant, resuspend the centrifuged cells in fresh RPMI-1640+10%FBS medium.
- 2) Inoculate the resuspended cells into a 96-well cell culture plate with a white wall and transparent bottom, 100ul/well of cell suspension, inoculate two culture plates, and place in a 37°C cell incubator for overnight culture.
- 3) The next day, take one of the 96-well plates inoculated with cells, add 100ul/well cell titer glo detection reagent and leave it for 60 minutes, read the value, and define it as G0 data.
- 4) Take another parallel plate, and take 11.1 ul/well of compound from



the previously diluted 10*compound concentration dilution plate to this plate , And set up another DMSO control well, continue to incubate in a 37°C cell incubator for 72 hours.

5) Take out the 96-well plate that has been treated with the compound for 72 hours from the incubator, add 100ul/well cell titer glo detection reagent and leave it for 60 minutes, read the value, and define it as G3 data.

6) Calculate the cell proliferation rate corresponding to each well according to the following formula %Proliferation = (Test compound well G3-G0 average value) / (DMSO control well G3 average value-G0 average value)*100.

7) According to the corresponding proliferation rate and concentration of each gradient concentration hole, use Prism Graphpad5.0 software to fit the gradient curve of cell proliferation, and calculate the GI50 of the compound (GI50 is defined as the corresponding compound concentration when the cell proliferation rate is 50% , The fitting formula in the software is as follows: $Y=Bottom + (Top-Bottom)/(1+10^{((LogIC50-X)*HillSlope)})$).

VIII. References

1. Congenital insensitivity to pain with anhidrosis (CIPA): effect of TRKA (NTRK1) missense mutations on autophosphorylation of the receptor tyrosine kinase for nerve growth factor. (PMID: 11159935) Mardy S ... Indo Y Human molecular genetics 2001
2. Congenital insensitivity to pain with anhidrosis: novel mutations in the TRKA (NTRK1) gene encoding a high-affinity receptor for nerve growth factor. (PMID: 10330344) Mardy S ... Indo Y American journal



of human genetics 1999

3. Mutation analysis reveals novel sequence variants in NTRK1 in sporadic human medullary thyroid carcinoma. (PMID: 10443680) Gimm O ... Eng C The Journal of clinical endocrinology and metabolism 1999 3

4. c.1810C>T polymorphism of NTRK1 gene is associated with reduced survival in neuroblastoma patients. (PMID: 20003389) Lipska BS ...

Limon J BMC cancer 2009

5. SNPs in neurotrophin system genes and Alzheimer's disease in an Italian population. (PMID: 18780967) Cozza A ... Pellegrini S Journal of Alzheimer's disease : JAD 2008

